


Application consideration


RCB0666 : HSNB

Please keep the following points in mind :

- Please do not change medium on the next day of the resuscitation (thawing and culture), and keep observing the cells for up to one week without medium change.
- The cells proliferate very slowly.
- Please be careful not to break up cell clumps too small when you passage them. Single cell culture is impossible for this cell line.

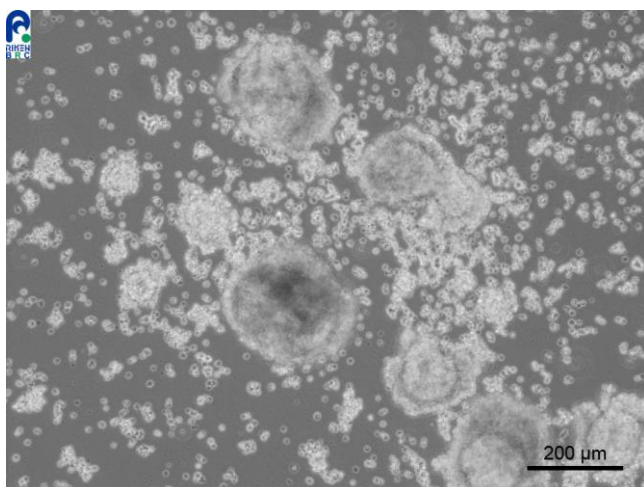

When cell clumps become too large, they start to turn dark color and result in cell death. Therefore, please passage them by breaking up the clumps before such condition.

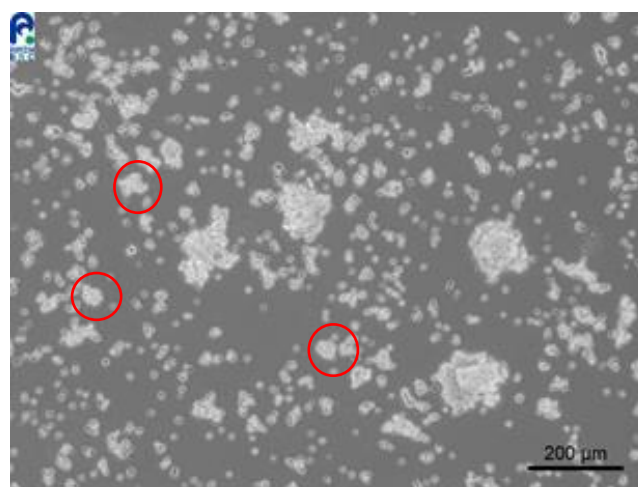

When you passage the cells, please break up cell clumps by pipetting to optimal sizes as shown in red circles.

【Technical support】 e-mail : cellqa.brc@riken.jp